


Eye Massager

Manual

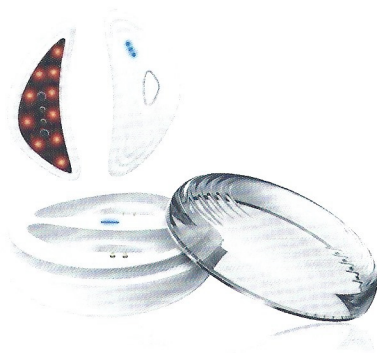

Please make sure to read this instruction manual carefully before use.

Preparation Before Use

- 1.Clean your face (remove makeup thoroughly), especially the eye bag area. Then apply a warm towel (ideally at 45°C) to the eye bag area for approximately 3 minutes.
- 2.Press and hold the main unit button to check if it enters the working state. If not, place the product in the charging case and charge it using the included USB data cable.


Product Description

Eye Beauty Instrument, This is an eye massage and care product that combines EMS microcurrent vibration massage with red light therapy. Through deep penetration via microcurrent vibration massage, it can promote the absorption of serums and provide eye care. It effectively relieves and improves eye concerns such as fine lines around the eyes, dark circles, sagging skin, and eye bags, while also soothing eye fatigue.

Precautions

- 1.If the gel eye patch loses moisture due to dry weather, you can use an appropriate amount of water spray to soften the internal molecular structure of the gel, restore its elasticity, and achieve a better massage effect.
- 2.A slight tingling sensation (like an electric touch) during use is normal. The intensity depends on skin moisture—higher moisture levels will result in a stronger sensation.
- 3.This product is not waterproof; do not clean it in water. If you need to clean it, ensure the device is turned off first, then wipe it with a damp tissue or towel.
- 4.It is recommended to use the product 3-4 times a week, 10 minutes each time. The specific usage duration may vary from person to person.
- 5.Store the product in a cool and well-ventilated place, away from direct strong sunlight and corrosive gases.
- 6.If the product is stored for more than 3 months, open the package to check it, recharge the device, and then use it.
- 7.Do not use the product if there is a wound on the eye area or if you have undergone eye surgery within one year.
- 8.Do not use this product in environments where the room temperature exceeds 35°C or drops below -5°C.
- 9.If there is any problem with the product, do not disassemble it by yourself. Please return it to the manufacturer for repair.

Product Introduction

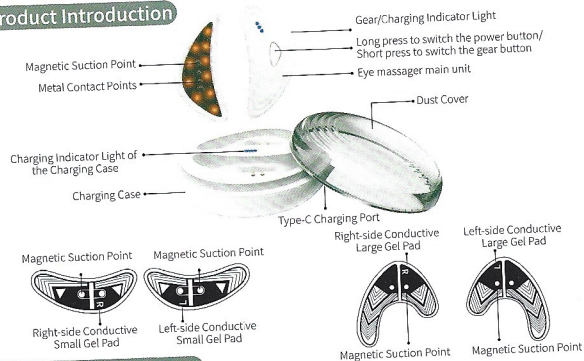


Product Specifications

Product Model: ES-008	Charging Voltage: 5V/1A	Power: $\approx 2W$
Gear: 3 Gears	Function: Red Light + EMS	Product Weight: 85g
Battery Capacity of the Eye Massager Main Unit: 27mAh	Battery Life: 120min	Charging Time: 180min
Battery Capacity of the Charging Case: 500mAh	Product Size: $\Phi 73 \times 36mm$	Package Size: $108 \times 108 \times 64mm$

Operating Mode

Gear 1: Red Light + Weak EMS	Gear 2: Red Light + Medium EMS	Gear 3: Red Light + Strong EMS
-------------------------------------	---------------------------------------	---------------------------------------

Operating Steps

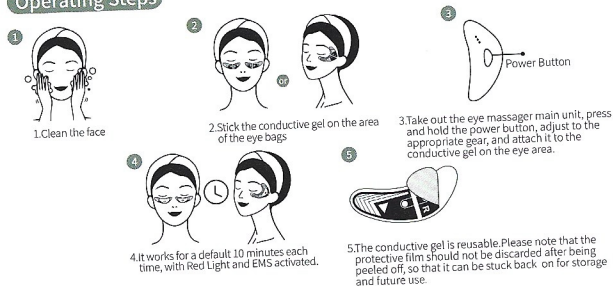

Frequently Asked Questions

Abnormal Phenomena	Reason	Solutions
Shuts down during use	Low battery level of the charging case	Please charge the product
Unable to charge	Incorrect voltage or current of the charger	Replace with a charger that meets the requirements in the manual for charging
No display during charging	Extremely low battery level	Please charge the product continuously
Long charging time	Extremely low battery level	Please charge the product continuously
No response during charging	Malfunction of the charger or charging cable	Replace the charger or charging cable

Note: To prevent malfunctions when the product is charging, do not use a 3-in-1 charging cable to charge this product together with other products at the same time.



When discarding this product, be sure to follow local regulations for disposal.